






Automated Specimen Processing, Robotic Workup and Digital Reporting: The Future of Microbiology

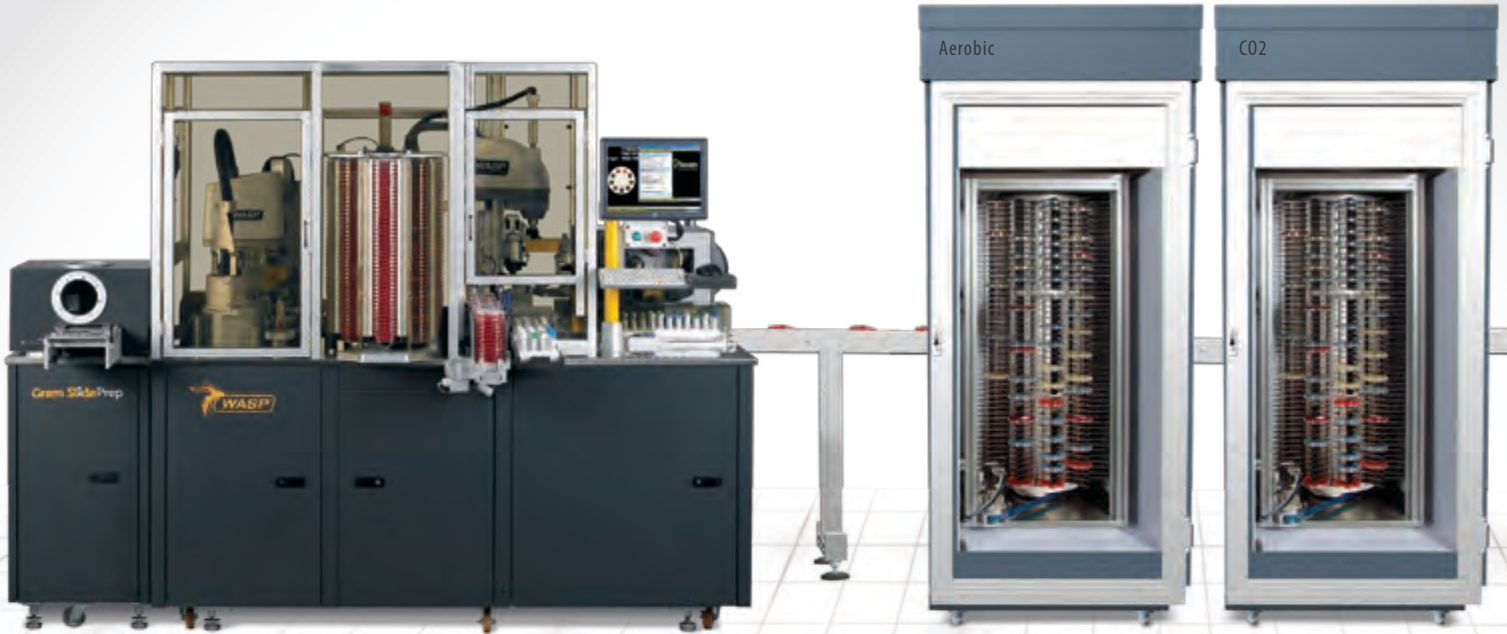

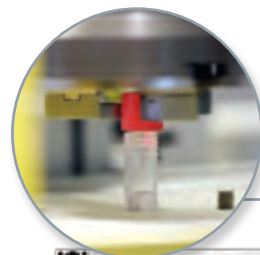

Automatically opens and closes any size specimen container

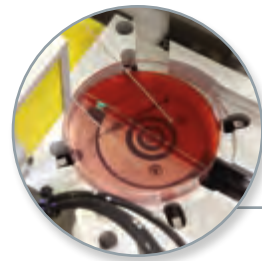

High throughput specimen management, planting and streaking

WASPLab™ is a barcode driven and conveyor connected specimen processing system utilizing robotic plate management and image analysis to automate specimen workup in Microbiology. Load and Walk Away, WASP® manages urine, swabs, sputum and feces directly, with no need for pre-processing.



Gram Slide Preparation

Automatic broth inoculation and Kirby Bauer disc dispensing



Universal Container Management System handles any container type individually.



ZEBRA Printer

Continuous feeding of any sample type
No Batching



Automatic loop and tool change station means longer walk away time and less user intervention

WASP® receives specimens into the LIS rather than manually for exact plating time

Nine silo carousel holds any manufacturer's plated media. Holds up to 378 plates



Smart Scan Technology reads barcode label regardless of position



User friendly touch screen monitor

Converging technology gives WASPLab™ the flexibility to merge seamlessly with new innovations.

Image Analysis

Analysis of digital culture plates can be done remotely. Robotic workup and digital reporting reduces turnaround time.




Robotic Incubator

Every plate has unique location for rapid access and homogeneous environmental conditions.


Automatic inversion of the plate prior to incubation.

Image Acquisition

scans plates at zero hour and at user defined intervals. Software can be set to compare before and after images and dispose of no growth plates.

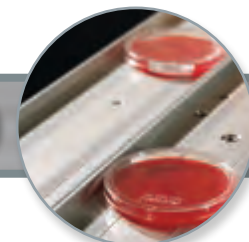

Conveyor fits any lab architecture and configuration

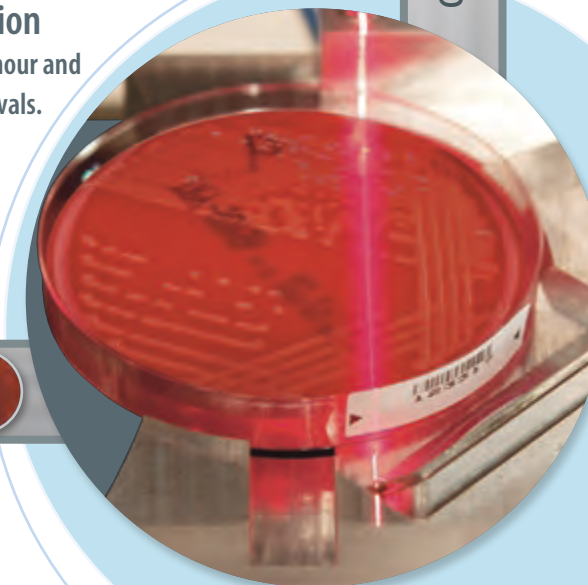

WASPLab™: The Future of Microbiology



innovation in preanalytics

Copan Italia, S.p.A.

Via F. Perotti 10, Brescia 25125, Italy

Tel: +39 030 268 7211 • Fax: +39 030 268 7250

e-mail: customer-care@copanitalia.com

website: copaninnovation.com

wasplab.com

Authorized Distributor